

● WAMIX-COMFORT IV, WAMIX-CLASSIC IV AND WAMIX-INJECTOR IV

It's all in the mixing!



Wamix-Comfort IV
with accessory
table stand

Wamix-Classic IV
with accessory
table stand

Wamix-Injector IV
with accessory
table stand

Technical information

Powerful, compact and innovative: The Wamix-mixers are an easy way of entering the world of vacuum mixing technology. The maintenance-free stirrer motor guarantees perfect, bubble-free mixing results with all kinds of plaster, embedding materials and silicone. The Wamix-Injector stands out with its powerful, oil-free and highly efficient vacuum generator, and the Wamix-Comfort / Wamix-Classic with its proven, maintenance-free motor pump. You could start stirring all machines by just pressing a button, the mixing time is adjustable up to 240 seconds or 999 seconds. Above that the Wamix-Comfort is equipped with a storage of up to 15 different programs with additional possibilities to store for example pre- and post-vacuum as well as other things.

- compact vacuum stirrers with powerful, maintenance-free motors
- perfect mixing of all kinds of plaster, embedding materials and silicone
- press-button start; automatic mixing when beaker is applied
- audible signal at programme end, time display flashes
- mixing time adjustable up to 240 seconds or 999 seconds
- beaker vented only after ON/OFF switch is activated
- separate motor and pump drive/vacuum generator drive, hence the agitator can be switched on and off at any time without interruption of vacuum

● Technical data and accessories

Wamix-Comfort IV	Item no.: 171977
Wamix-Classic IV	Item no.: 171979*
Wamix-Injector IV	Item no.: 171978**
Voltage*	220–240 V / 50/60 Hz
Power consumption	2.9 A (Si 4,0 A), 2.4 A (Si 4,0 A)*, 1.8 A (Si 4,0 A)**
Output	315 W, 220 W*, 120 W**
W x H x D incl. wall mount	205 x 370 x 245 mm
Weight	9.9 kg, 9.5 kg*, 8.3 kg**
Sound level	≤ 70 dB (A)
Speed	200 - 600 min ⁻¹ , 200 - 600 min ⁻¹ *, 380 min ⁻¹ **
Vacuum	approx. 900 mbar, approx. 900 mbar*, approx. 950 mbar**
Mixing time	up to 240 sec., 999*/**sec. adjustable
Safety class	IP 22
compressed air connection	min. 5 bar**

*Other voltages on request

Included in delivery	Item no.:
Beaker, medium cpl. 350 ml	171283
Mounting set	174358
VC hose Ø 8 x 6 textile (blue) 2 m**	320026
Accessories	Item no.:
Table stand (W 219 x H 467 x D 308 mm)	171777
Beaker, extra small cpl. 75 ml	171295
Beaker, small cpl. 200 ml	171230
Beaker, medium cpl. 350 ml	171283
Beaker, large cpl. 600 ml	171237
Beaker, extra large cpl. 1000 ml	171269
PVC hose Ø 8x6 textile (blue) per meter**	320026
Vacuum washing beaker 600 ml	171275


Technical changes reserved



Wamix-Comfort IV with wall mount



Wamix-Classic IV and Wamix-Injector IV with wall mount

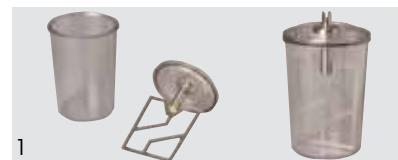

1. 1000 ml, 2. 600 + 350 ml, 3. 75 ml + 200 ml

Major features at a glance

- compact vacuum stirrers for low-cost entry to vacuum technology
- bubble-free mixing of plasters, embedding materials and silicone
- user-friendly operation, push-button start
- low-maintenance
- simple cleaning thanks to stay-clean surfaces
- 350 ml beaker and wall bracket supplied, table stand and other beaker sizes available as extras

www.wassermann-dental.com